



# 6050

## Calorimetro Isoperibolico

---

### **Descrizione**

Il nuovo calorimetro Parr 6050 è uno strumento che utilizza una tecnica di misura compensata in continuo. Questo modello è stato progettato per essere gestito da un Personal Computer.

### **Caratteristiche**


- Controllo da PC mediante collegamento diretto a un PC Windows con cavo USB standard. Tutte le funzioni e le prove sono accessibili con un semplice clic del mouse utilizzando il software dedicato;
- Riempimento automatico dell'ossigeno: è caricata la giusta quantità di ossigeno per ogni test;
- Lo strumento lavora a temperatura ambiente;
- Il tempo per ogni test è di circa 8 minuti in modalità dinamica e 12 minuti in modalità equilibrio.

Parr Instrument Company



# 6050 Compensated Jacket Calorimeter

**Designing  
and Building  
High Precision  
Calorimeters  
for Over  
100 Years**

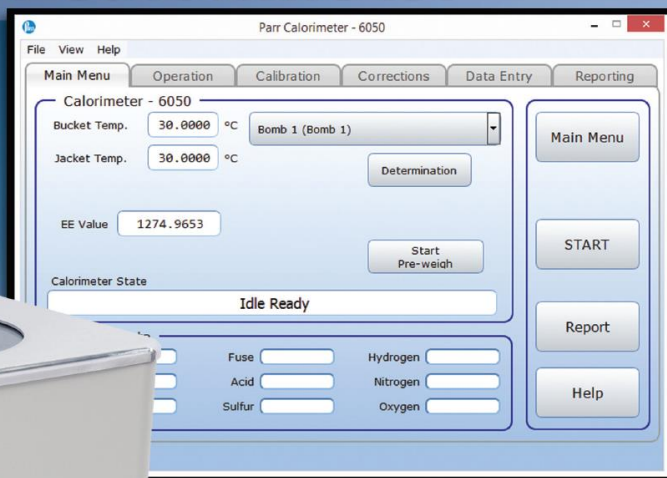

Sample screenshot from a PC running the Parr 6050 Calorimeter Software.

**T**he new Model 6050 Compensated Jacket Calorimeter is the first calorimeter designed by Parr to be controlled by your own personal computer, helping to make it Parr's most affordable calorimeter to date.

**Features of the new 6050 Calorimeter include:**

- **PC Control:** The 6050 is attached directly to a Windows PC using a standard USB cable. All functions and tests are accessible with the click of a mouse using Parr's 6050 Software Package, installed on your own PC
- **Automatic Oxygen Filling:** The 6050 fills the combustion vessel with the correct amount of oxygen for each test automatically – removing one manual step from the standard combustion vessel preparation routine
- Operates at room temperature
- Operator time per test: 6 minutes

**Let Us Build One for You.**